

# Specification Sheet

Liquid Chromatograph Mass Spectrometer

## LCMS-8060

The LCMS-8060 delivers the highest sensitivity and fastest analysis speed of any LCMS on the market today. A newly developed UF-Qarray boosts ion intensity but suppresses noise. By improving the ion sampling device, the ion guide, and vacuum efficiency, Shimadzu has achieved an unprecedented sensitivity in quantitative analysis by LC/MS/MS. Shimadzu's proprietary ultra-fast technologies (UF Technologies) allow acquisition of high-speed MRM transitions, ultra-fast polarity switching, and ultra-fast scanning, all with high data quality. In conjunction with a robust design and easy maintenance, the LCMS-8060 offers the ideal solution for better data quality and higher throughput for a wide array of challenging and routine applications.



### Instrument

Model	LCMS-8060
Mass range	$m/z$ 2 to 2,000
Sensitivity	ESI positive: 1 pg reserpine, S/N > 750,000:1 (RMS) ESI negative: 1 pg chloramphenicol, S/N > 750,000:1 (RMS)
Resolution	$R < 0.7 \mu$ FWHM
Mass stability	0.05 $\mu$ /12 hr
Cross-talk	< 0.003 %
Minimum pause time	1 msec
Minimum dwell time	0.8 msec
Scan speed	Max 30,000 $\mu$ /sec (in all modes of scanning) (0.1 $\mu$ step: 300,000 data points/sec)
Polarity switching time	5 msec
Interface	ESI (Standard), APCI (Optional), DUIS (Optional)
Applicable LC flow rate	ESI 1 $\mu$ L/min to 2 mL/min
MRM transition speed	Max 555 channels/sec
DL maximum temperature	300 °C
Block heater maximum temperature	ESI/DUIS: 500 °C APCI: 300 °C
Interface maximum settable temperature <sup>Note)</sup>	ESI/DUIS: 400 °C APCI: 500 °C

Analysis mode	Q1 Scan/SIM Q3 Scan/SIM MRM Precursor ion scan Product ion scan Neutral loss scan
---------------	--

### Mass Analyzers and Detector

Mass analyzers	Q1 & Q3 are molybdenum hyperbolic mass filters with pre-rods; Q1 includes post-rods
Collision cell	Tapered multipole type ultra-high-speed collision cell (UFSweeper III collision cell)
Detector	Secondary electron multiplier with off-axis conversion dynode
Ion optics	Q-array focus optics operating in Field-Flow mode, multipole transfer optics
Digital detection system	Operates in pulse counting mode for fastest operation
Detection mode	Ultra-fast positive/negative ion switching
Dynamic range (pulse counting)	$2 \times 10^7$ cps
Vacuum system	Rotary pump: 1 unit Vacuum pumping speed: 65 $m^3$ /hour Triple-inlet turbo molecular pump: 1 unit 40 L/sec, 260 L/sec, 210 L/sec


Note) Actual maximum ESI/DUIS temperature when set at 400 °C on LabSolutions LCMS is greater than 700 °C.  
Actual maximum APCI temperature when set at 500 °C on LabSolutions LCMS is greater than 600 °C.

## Software


Workstation	LabSolutions LCMS
Operating system	Windows 7 32 Bit, 64 Bit
Instrument control	Prominence and Nexera series
MS acquisition mode	Scan (Max. 1,000 events), SIM (Max. 1,000 events × 32 channels)
MS/MS acquisition mode	MRM (Max. 1,000 events × 32 channels) Product ion scan Precursor ion scan Neutral loss scan
Auto-tuning	Possible to optimize sensitivity, resolution, and mass calibration in both positive and negative ionization mode

## Installation Example

[Plan View]

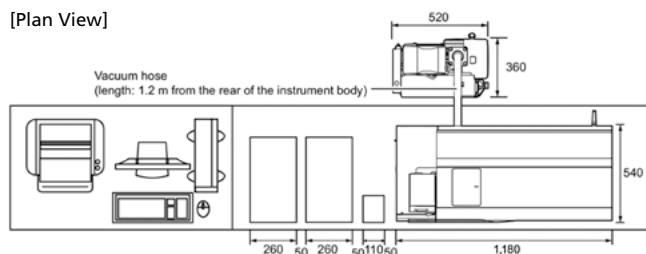

[Front View]

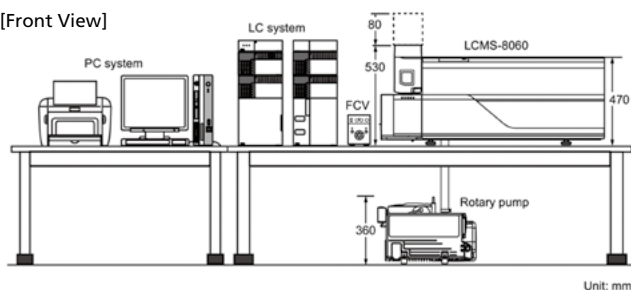

## Installation Conditions

Temperature	18 to 28 °C
Humidity	20 to 70 % (Non-condensing)
Size	1,180 mm (W) × 540 mm (D) × 610 mm (H)
Weight	140 kg
Power supply	MS unit: AC 230 V 15 A (50/60 Hz) Single-phase Rotary pump: AC 200–240 V 12 A (50/60 Hz) Single-phase
Gas requirements	Nitrogen gas: Maximum 24.4 L/min, Purity greater than 97 % Argon: Purity greater than 99.99% as CID gas Dry air: Maximum 20 L/min, oil/water-free Total nitrogen plus air: 25L/min maximum

The above are not standard installation specification. All LCMS-8060 instruments will be installed and tested in accordance with standard performance tests as detailed in the Shimadzu document ZEAH-0486, Shimadzu High-Performance Liquid Chromatograph Mass Spectrometer LCMS-8060 Installation Standard.



Shimadzu Corporation

[www.shimadzu.com/an/](http://www.shimadzu.com/an/)

**For Research Use Only. Not for use in diagnostic procedures.**

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®". Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®". Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.

© Shimadzu Corporation, 2017